
Infrared Ear Thermometer

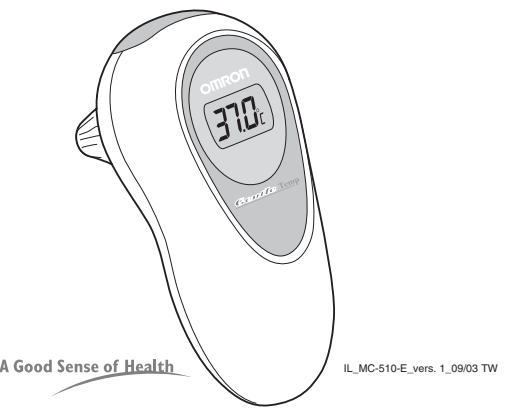Important instructions for accurate temperature measurement
Wichtige Anweisungen für eine akkurate Temperaturmessung
Belangrijke aanwijzingen voor nauwkeurige temperatuurmeting
Instructions importantes pour une mesure exacte de la température
Istruzioni importanti per una misurazione accurata della temperatura
Instrucciones importantes para una precisa medición de la temperatura

Introduction
This leaflet contains extra information to achieve accurate temperature measurements. Read the OMRON Instruction Manual and follow the instructions on this leaflet carefully, before you use the OMRON infrared Ear Thermometer.

- General advice**
- 1 Familiarize yourself with the ear thermometer by practicing on yourself and family members. By doing this you improve your handling technique and tell when a reading is higher than normal.
 - 2 When using the thermometer for the first time, or when an accurate reading is essential, take three measurements and note the highest value. Always use the same ear when comparing readings.
 - 3 When applied on small children cover the ear opening gently with the tip of the probe while taking the measurement.
 - 4 Do not rely on temperature measurements alone. If you have health concerns, seek medical advice.

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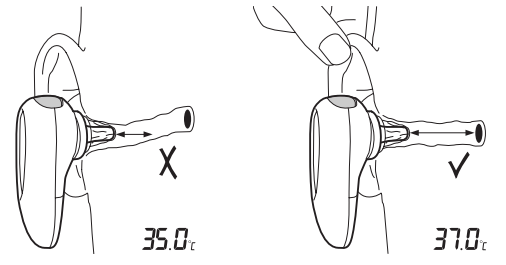
The OMRON Infrared Ear Thermometer may give low readings when not correctly placed in the ear. The thermometer detects heat from the eardrum, but also responds to heat from the ear canal. The temperature of the ear canal could be lower (approximately 2°C) than the eardrum.

There are two modes to measure temperature of the eardrum. Mode A takes 1-3 seconds to measure, mode B takes 10 seconds.

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